

# PROTECTION THAT WORKS FOR YOU AND YOUR PATIENTS

PPE serves a dual purpose. Originally developed to protect healthcare workers, PPE, when used properly, can also represent the first line of defense against contact transfer to patients of pathogens such as MRSA and C. diff.

Our AAMI Level 2 Over-the-Head Isolation Gown with Thumb Hooks meets or exceeds the AAMI PB70 requirements for AAMI Level 2 fluid protection.

When you don fresh PPE before entering each patient's room, you help reduce the risk of contamination and transmission of infectious organisms that can lead to HAIs.



### **AAMI LEVEL 2 OVER-THE-HEAD ISOLATION GOWN WITH THUMB HOOKS**

- Meets AAMI Level 2 requirements
- Medium-weight Tri-layer fabric
- Ties at side for easier gown closure.

### **RECOMMENDED FOR:**

- Suturing, blood draw, inserting I.V. lines, specimen handling, drawing arterial blood
- ICU, Hyperbaric and Dialysis Units, Labs and Pathology, Nursery

## AAMI PB70 Guidelines: A Reliable Benchmark for Protection

	LEVEL 1	LEVEL 2	LEVEL 3
AAMI Level / Performance Requirements at 4% AQL	Spray Impact Penetration $\leq$ 4.5g	Spray Impact Penetration < 1.0g Hydrostatic Pressure > 20cm	Spray Impact Penetration < 1.0g Hydrostatic Pressure > 50cm
Gown Characteristics (Weight)	Light-weight SMS fabric	Medium-weight SMS fabric	Heavy-weight SMS fabric

To meet AAMI Guidelines, protective apparel must achieve minimum performance standards for strength, barrier performance and fluid resistance. Tested under controlled conditions, AAMI-rated gowns must demonstrate a 4% AQL level, which provides a high level of assurance that users are protected for the anticipated level of fluid contact. Each garment must be fully protective in the front, back and along the seams, and be clearly labeled to indicate the level of protection provided, from AAMI Level 1 (lowest) to AAMI Level 4 (highest). The Halyard Tri-layer AAMI Level 2 Over-the-Head Isolation Gown with Thumb Hooks meets all AAMI specifications for Level 2 protection.

### HALYARD\* TRI-LAYER AAMI LEVEL 2 OVER-THE-HEAD ISOLATION GOWN WITH THUMB HOOKS


Code	Description	Cuff End	Size	Gowns/Case	Bags/Case
43146	Thumb Hooks, Yellow, Closed Back	Thumb Hooks	Large	100 Case	10
43147	Thumb Hooks, Yellow, Closed Back	Thumb Hooks	X-Large	100 Case	10
44715	Thumb Hooks, Yellow, Closed Back	Thumb Hooks & Elastic	Large	100 Case	10
44717	Thumb Hooks, Yellow, Closed Back	Thumb Hooks & Elastic	X-Large	100 Case	10
44719	Thumb Hooks, Yellow, Closed Back	Thumb Hooks & Elastic	XX-Large	100 Case	10

### HALYARD\* TRI-LAYER OVER-THE-HEAD ISOLATION GOWN WITH THUMB HOOKS

Code	Description	Size	Fabric Weight	Gowns/Case	Bags/Case
29549	Thumb Hooks, Yellow, Perforated Open Back	Universal	Medium	100 Case	10
29516	Thumb Hooks, Yellow, Perforated Open Back	XX-Large	Medium	100 Case	10
32574	Thumb Hooks, Yellow, Perforated Open Back, 11" x 5.5" Fold	Universal	Medium	120 Case	12
32572	Thumb Hooks, Yellow, Perforated Open Back, 11" x 5.5" Fold	XX-Large	Medium	120 Case	12

### HALYARD\* BASICS\* TRI-LAYER OVER-THE-HEAD ISOLATION GOWN WITH THUMB HOOKS

Code	Description	Size	Fabric Weight	Gowns/Case	Bags/Case
43148	Thumb Hooks, Yellow, Perforated Open Back	Universal	Light	150 Case	6
43149	Thumb Hooks, Yellow, Perforated Open Back	X-Large	Light	150 Case	6

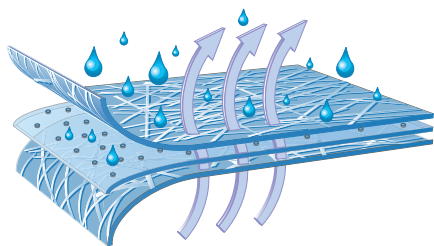

The HALYARD\* Over-the-Head Isolation Gown with Thumb Hooks is Made with Protective Tri-Layer Fabric

Outer layers—Spunbond polypropylene material that is cloth-like, yet tough and abrasion-resistant

Middle layer—Intricate meltblown polypropylene micro fibers that filter out many water-based liquids



KNOWLEDGE NETWORK\* Accredited Education  
In-Service Training and Technical Support  
Credentialed Sales Representatives  
Tools & Best Practices  
Clinical Research  
Commitment to Excellence

[www.halyardhealth.com](http://www.halyardhealth.com)

