

HALYARD*
Ostomy Care



YOUR PARTNER IN SKIN AND OSTOMY CARE

BECAUSE SKIN TAKES CARE





BECAUSE SKIN TAKES CARE

The skin plays a crucial role in overall health by acting as a protective barrier to bacteria, chemicals, external forces and excessive moisture. One of the leading causes of skin breakdown is moisture-related skin damage caused by extended exposure to bodily fluids, putting patients at risk for inflammation or even infection.

Our line of ostomy care products provides essential protection for patients living with a stoma, and support for those who care for them. Specialized barrier seals, pouch deodorants and wipes protect delicate skin from damage caused by moisture and adhesive stripping. Intact, healthy skin means better fitting ostomy products and less chance of leakage.

Protect with confidence when you protect with O&M Halyard.



Why Use an Ostomy Care Accessory?

Ostomy accessories work in tandem with your pouching system to provide various benefits. Here's why you might want to consider incorporating them into your routine:

- **Enhance Pouch Performance:** Accessories can extend wear time, safeguard your skin, and simplify adhesive removal.
- **Address Common Issues:** They can help resolve problems like leakage and skin irritation.
- **Improve Quality of Life:** By eliminating odors, accessories can significantly enhance your daily comfort.

Every individual and stoma is unique. Barrier rings play a vital role in achieving optimal results for patients:

- **Customizable Fit:** These rings can be stretched, shaped, or cut into pieces to fill uneven areas.
- **Leakage Prevention:** They help prevent stoma output from seeping under the pouch seal, reducing leakage and skin irritation.
- **Flexible Use:** The rings can be cut or stacked to better fit your pouching system.
- **Gentle on Skin:** Being alcohol-free, they cause no stinging when applied to broken or irritated skin.

These features ensure that ostomy accessories not only address practical concerns but also contribute to a more comfortable and confident lifestyle.



HALYARD* Barrier Seals — Rings, Arcs and Sheets

Rings, Arcs and Sheets designed to protect the peristomal skin and aid in securement of wafer

The skin around a stoma is at risk for skin breakdown caused by effluent from the stoma, and/or repeated removal of ostomy wafer which can strip the skin. Ostomy patients with unique stoma shapes or irregular peristomal skin which can be 'filled' with these flexible, moldable barrier rings.

Ostomy Barrier Seal Rings can be shaped to fit the stoma and applied either to the peristomal skin or to the wafer before application to the stoma site. Precut starter holes can be cut and stretched to fit the shape of the stoma.

Ostomy Rings are available in three styles: Conforming (stretchable), Intact (low residue/no-meltdown formula) and Dual Layer (two sides have variable absorbency).

Ostomy Barrier Arcs can be used under or over the ostomy wafer to help improve the fit around the stoma or to window frame the edges to aid in a more leak-proof edge of wafer. The barrier arcs can help protect the skin from adhesive damage, conforms to fit unique shape of stoma or skin surface and provides a water-resistant cover to address leakage concerns.

Ostomy patients with unique stoma shapes or irregular peristomal skin which can be 'filled' with these flexible, moldable arcs. Available: 1/3 circle, 1/3 circle extra wide, 1/2 circle and 4x4 sheet.

Key Features of Barrier Rings, Arcs and Sheets

- Skin-friendly adhesive barriers can be used under wafer or over wafer to secure edge.
- Absorbent, protects skin from moisture or adhesive stripping and helps extend wafer wear-time
- Thin, comfortable, easy removal adhesives, little to no residue



HALYARD* Ostomy Powder

Ostomy powders are used to dry a raw, weepy area on the skin. The powder absorbs the moisture from the skin allowing the bag to adhere properly, especially if the skin is very moist.

Application:

1. Make sure the skin surrounding the stoma is clean, dry and free from oils, creams, and/or soap residues. Untwist the cap to open powder.
2. Apply directly on skin: Dust a light coating of powder onto the affected area, brush off excess.
 - a. For creating Crusting Technique: Dust the powder onto skin. Use a skin barrier wipe or spray over the top to create a crust to protect traumatized skin.
3. Allow to completely dry before adhering appliance. Available 1 oz. and 5 oz. Puff bottle for convenient application.


HALYARD* Ostomy Deodorant

Deodorant Drops to Eliminate Odor in Ostomy Pouches

Ostomy deodorants are not body deodorants but intended to help manage odor in the ostomy pouch. Odor eliminator drops are placed inside the ostomy pouch and help eliminate odor when pouch is emptied. Use with each pouch change or after emptying. Available in clear, unscented with lubricant to aid in emptying pouch contents. Also available in blue, unscented.

Key Features:

- Available clear or blue color liquid
- Unscented, Helps Eliminate Odor
- Available with Lubricant to help empty pouch contents

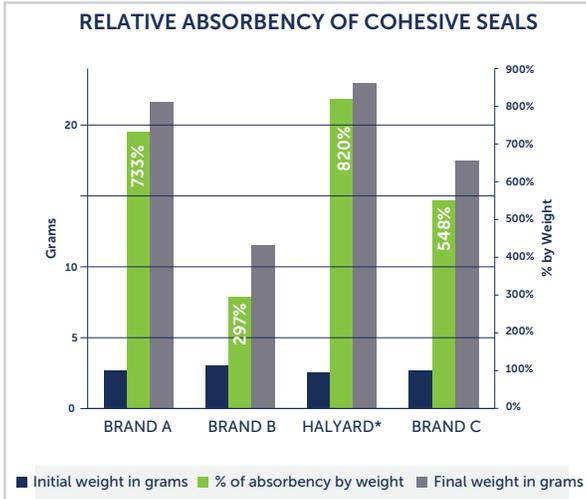


Source: HALYARD* data on file.

Absorbency Comparison: Leading brands of comparable Skin Barrier Seals were tested to show absorbency performance.

Results: HALYARD* Barrier Seals outperformed the competition in absorbency performance of total grams and % by weight.


Relative Tack of Cohesive Seals: Leading brands of comparable Skin Barriers were tested to show relative tack.

Results: HALYARD* Skin Barriers outperformed the competition in tack of cohesive seals.

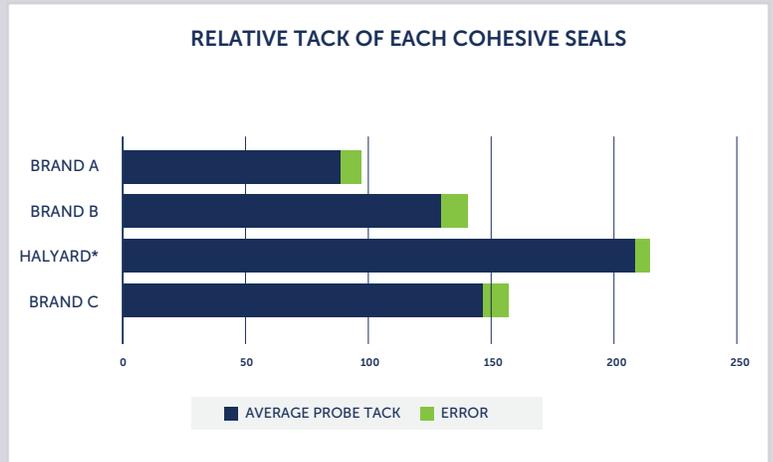

Source: HALYARD* data on file.

KEY ATTRIBUTE OF MOLDABLE SEALS TURTLE-NECKING AROUND STOMA

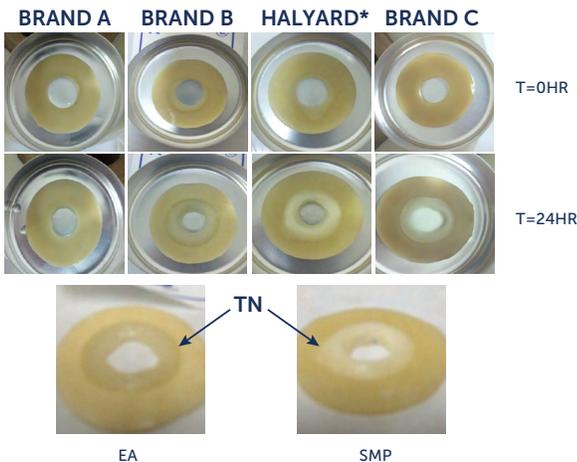

Conformability and Turtle-Neck Feature: Leading brands of comparable Skin Barrier Seals were tested to show conformability and turtle-necking feature.

Results: HALYARD* Barrier Seals have comparable conformability and turtle-neck performance to other brands of seals.

Source: HALYARD* data on file.

Product Code	Description	Size	Qty / Sales UM	Base UoM	Sales UoM	Qty / Inner	Inner Pkg
Barrier Rings							
55601	HALYARD* Barrier Ring Conforming	2in. OD	160	EA	CA	10	BX
55602	HALYARD* Barrier Ring Conforming	2in. OD	400	EA	CA		
55616	HALYARD* Barrier Ring Conforming	2in. OD 4.2mm Thick	400	EA	CA		
55619	HALYARD* Barrier Ring Conforming	4in. OD	160	EA	CA	10	BX
55671	HALYARD* Barrier Ring Intact Seal	2in. OD	160	EA	CA	10	BX
55678	HALYARD* Barrier Ring Dual Layer	2in. OD	160	EA	CA	10	BX
Barrier Arcs & Sheets							
55636	HALYARD* Barrier Crescent 30/Pk	1/3 Circle Crescent	1,500	EA	CA	30	BX
55661	HALYARD* Barrier C-Shape Strips 20/Pk	1/2 Circle Crescent	1,000	EA	CA	20	BX
55666	HALYARD* Barrier Crescent Extra Wide 30/Pk	1/3 Circle Crescent Extra Wide	900	EA	CA	30	BX
55684	HALYARD* Barrier Sheet	4in. x 4in.	50	EA	CA	5	BX
Ostomy Deodorant							
55644	HALYARD* Ostomy Deodorant Lubricate	.25 oz. Sachet	400	EA	CA		
55645	HALYARD* Ostomy Deodorant Lubricate	8 oz. Bottle	12	EA	CA		
55658	HALYARD* Ostomy Deodorant	8 oz. Bottle	12	EA	CA		
Barrier Powder							
55755	HALYARD* Barrier Powder Ostomy	5 oz. Puffer Bottle	40	EA	CA		
55754	HALYARD* Barrier Powder Ostomy	1 oz. Puffer Bottle	40	EA	CA		

www.halyardhealth.com



HALYARD* is an Owens & Minor brand

**Talk with your O&M Halyard
representative to get more details
and order these products.**

*Registered Trademark or Trademark of O&M Halyard or its affiliates. ©2025. All rights reserved.
008865 COPY-05406