

MediChoice®

# WASTE FLUID MANAGEMENT SYSTEM

Disposable Suction Canisters

Suction Tubing & Yankauers



NOT MADE WITH  
natural rubber latex

sterile EO

MEDIChoice®

# suction canisters


disposable

## Key Features

The MediChoice® Suction Canister System is single-use, designed for convenience by removing the need for a secondary waste liner. The canisters include an anti-foaming emulsion that prevents froth formation and aerosolization, reducing the risk of contaminant exposure while maximizing fluid capacity.

MediChoice® Canisters support best practice by integrating a non-mechanical shut-off valve that prevents back flow of aerosols and fluids into the vacuum system, and includes a 0.3 micron bacterial eliminating filter effective in removing 99.9% of pseudomonas.

- Lid attaches firmly to the canister, maximizing vacuum power and reduces risk of leakage
- Easy-to-read volume marks help with patient documentation
- Reflux protector on patient port reduces contamination exposure risk

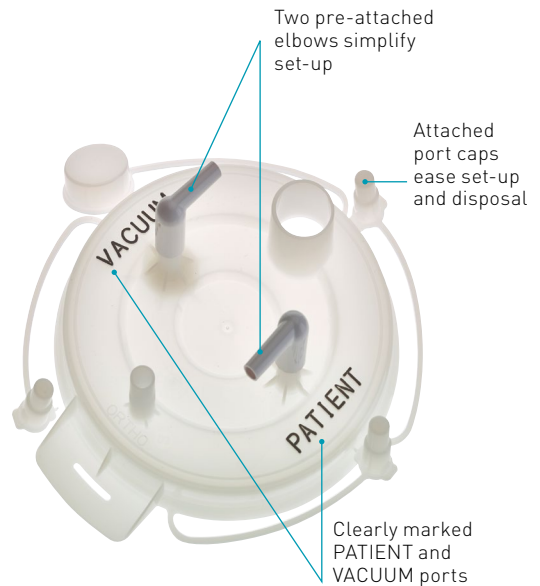

MediChoice® Disposable Suction Canisters			
SKU	Description	Size	Quantity
WFM0700D	700 mL Disposable Suction Canister with Lid and 2 Elbows	700 mL	32 EA/CA
WFM1200D	1200 mL Disposable Suction Canister with Lid and 2 Elbows	1200 mL	32 EA/CA
WFM2000D	2000mL Disposable Suction Canister with Lid and 2 Elbows	2000mL	24 EA/CA
WFM3300	3300 mL Disposable Suction Canister with Lid and 2 Elbows	3300 mL	12 EA/CA


# suction tubing

with connectors

## Key Features

Suction Tubing, with molded female connectors and male adapters, provides secure connections to support uninterrupted procedures. The striated PVC tubing with low coil memory is strong and collapse-resistant at high vacuum pressure yet flexible for ease of draping, promoting best clinical outcomes.

- Sterile and non-sterile tubing
- Inner diameter sizes: 3/16 in., 1/4 in., and 9/32 in.
- Tubing lengths: 6 ft., 10 ft., 12 ft. and 20 ft.

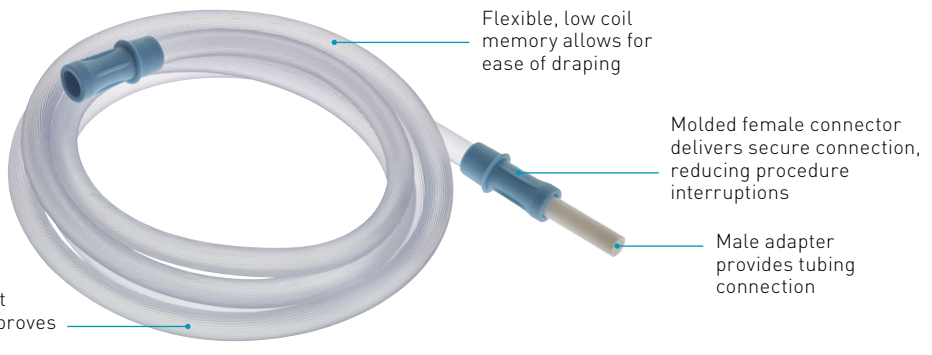


SKU	Description	Quantity
ST1001	Sterile, 6 Feet x 3/16 Inch, green connector	50 EA/CA
ST1002	Sterile, 10 Feet x 3/16 Inch, green connector	50 EA/CA
ST1003	Sterile, 12 Feet x 3/16 Inch, green connector	20 EA/CA
ST1004	Sterile, 20 Feet x 3/16 Inch, green connector	20 EA/CA
ST2001	Sterile, 6 Feet x 1/4 Inch, light blue connector	50 EA/CA
ST2002	Sterile, 10 Feet x 1/4 Inch, light blue connector	50 EA/CA
ST2003	Sterile, 12 Feet x 1/4 Inch, light blue connector	20 EA/CA
ST2004	Sterile, 20 Feet x 1/4 Inch, light blue connector	20 EA/CA
ST3001	Sterile, 6 Feet x 9/32 Inch, medium blue connector	50 EA/CA
ST3002	Sterile, 10 Feet x 9/32 Inch, medium blue connector	50 EA/CA
ST3003	Sterile, 12 Feet x 9/32 Inch, medium blue connector	20 EA/CA
ST3004	Sterile, 12 Feet x 9/32 Inch, medium blue connector	20 EA/CA

# suction yankauers

## Key Features

Suction Yankauers are designed with shatter-resistant materials. The slip-resistant, lightweight handle helps provide optimal control and maneuverability. A translucent pathway allows for fluid visualization and the smooth tips are designed with side eyes to minimize the chance of tissue grab and damage.

- Sterile vented and non-vented
- Rigid and flexible handles
- Regular and Bulbous tips

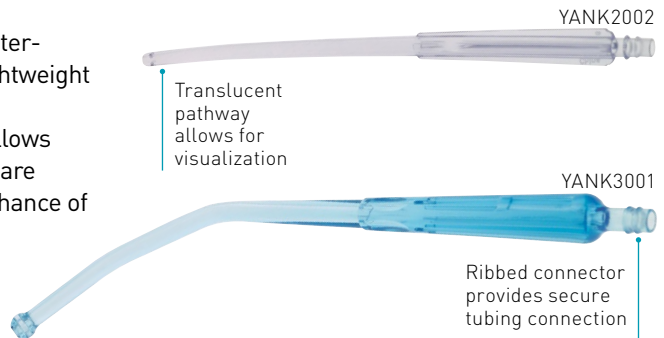

### MediChoice® Suction Yankauers

SKU	Description	Quantity
YANK1001	Suction Yankauer with Bulb Tip, Non-Vented, Flexible, Polyvinyl Chloride, Blue, Sterile	50 EA/CA
YANK1002	Suction Yankauer with Bulb Tip, Vented, Flexible, Polyvinyl Chloride, Blue, Sterile	50 EA/CA
YANK3001	Suction Yankauer with Bulb Tip, Non-Vented, Rigid, Polyvinyl Chloride, Purple, Sterile	50 EA/CA
YANK2001	Suction Yankauers with Open Tip, Non-Vented, Flexible, Polyvinyl Chloride, Blue, Sterile	50 EA/CA
YANK2002	Suction Yankauers with Open Tip, Vented, Flexible, Polyvinyl Chloride, Blue, Sterile	50 EA/CA

# MEDI CHOICE®

---

To order MediChoice Waste Fluid Management Systems or any other MediChoice product, please contact your local Owens & Minor representative at +1 (800) 488-8850, or visit us at [medichoice.com](http://medichoice.com)

---

**OWENS & MINOR, INC** (NYSE: OMI) is a global healthcare services company dedicated to *Connecting the World of Medical Products to the Point of Care<sup>SM</sup>* by providing vital supply chain services to healthcare providers and manufacturers of healthcare products. Owens & Minor provides logistics services across the spectrum of medical products from disposable medical supplies to devices and implants. With logistics platforms strategically located in the United States and Europe, Owens & Minor serves markets where three quarters of global healthcare spending occurs. Owens & Minor's customers span the healthcare market from independent hospitals to large integrated healthcare networks, as well as group purchasing organizations, healthcare products manufacturers, and the federal government. A FORTUNE 300 company, Owens & Minor is headquartered in Richmond, Virginia, and has annualized revenues exceeding \$9 billion. For more information about Owens & Minor, visit the company website at [www.owens-minor.com](http://www.owens-minor.com), follow [@Owens\\_Minor](https://twitter.com/Owens_Minor) on Twitter, and connect on LinkedIn at [www.linkedin.com/company/owens-&-minor](http://www.linkedin.com/company/owens-&-minor).

MediChoice is a registered trademark of Owens & Minor, Inc.

© 2017. All rights reserved.